

FINE MARKETING

15A, BHARAT INDUSTRIAL ESTATE, LAKE ROAD, BHANDUP(W), MUMBAI-400078.INDIA.

TEL: 022-25965796 / 95 FAX: 022- 25964859

EMAIL: siyer@bom5.vsnl.net.in / finemarketingindia@gmail.com

Website : www.finemarketingindia.com / www.materialtestingmachines.net

Cold setting Resin

When the specimen is of odd shape or size to hold it while polishing it becomes essential to mount it.
(One pack Consists Liquid 400 ml + powder 400 gm)

Cold setting resin




Die Set

For mounting the specimens with above cold setting
They are kept face down in the die set & cold setting
resin in 1:1 powder liquid proportion is poured in cavity.
Mounted specimens become ready in five minutes.
Mould diameter 25,32,38 & 50 mm. Thickness 15mm

Die set

Bakelite mounting press

For mounting the specimens with Bakelite powder
Specimen is kept face down in the Cavity & bakelite
Powder is Hot moulded around it with hydraulic
pressure. Mould dia 30 mm.
With Digital timer. Operates on 220 V A/C

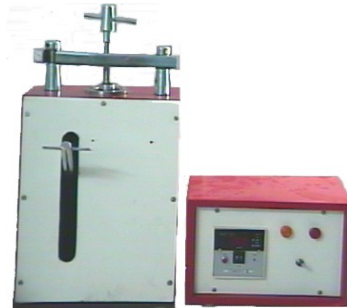

Bakelite powder (10 kgs/pack)